



HK-100 Veterinary Infusion Pump

Feature & Functions:

Compact design

The compact and light weight design saves space
Large & colorful LCD display

Added safety

Free flow protection

The built-in flow stop feature helps reduce the risk of accidental free flow when pump door is open.

Anti-bolus function

Diminishes the volume of unwanted Bolus injected to the patient after removal of the occlusion cause.

Anti-reverse rotation detection

Panel Lock

The panel lock feature helps prevent unauthorized changes of any instrument setting

Versatile functions

2000 history records

Real time display

RS232 interface (optional)

90° rotation pole clamp convenient for horizontal bar, vertical IV pole & platform use

User-friendly operation

Soft key design, easy to operate

Load directly the last infusion rate and volume limit

Adjustable buzzer volume

Top indicator light reflects the operation status as well as alarm case

Night mode to ensure good rest of patients

Removable pump body for easy cleaning



Free flow protection




Rs232 interface

2000 event logs and alarm information can be transferred to computer for review



90° rotatable pole clamp

Fixation on vertical pole or horizontal bar

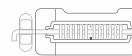

Removable pump body

convenient for cleaning

HK-100 Veterinary Infusion Pump

Specification

Flow rate range	0.1 ~ 1200ml/h Micro mode: 0.1-99.9 ml/hr in 0.1 ml increments
Volume to be infused (VTBI)	0 ~9999ml
Volume Infused (Σ)	0 ~36000ml
Infusion accuracy	$\pm 5\%$
Applicable infusion set	Support a wide variety brands of infusion set including major manufacturers & selected local brands
KVO rate	1 ~ 5 ml/h
Bolus rate	300 ~ 1200 ml/h
Purge	600ml/h
Air bubble detection	Ultrasound sensor detection with four levels of air in line alarm adjustable
Occlusion pressure	3 adjustable occlusion pressure settings: low, middle and high
Alarms	Visual and audible alarms:Door open,Air-In-Line, Occlusion,Infusion completion, infusion near over, Start reminder function, Low Battery,Battery depleted,malfunction etc.
Computer Interface	Rs232 (optional)
History records	2000 history records
Power supply	AC:100~240V,50/60Hz DC:12V $\pm 1V$
Battery	Rechargeable lithium polymer battery,7.4V,1900mAh Can operate for more than 3 hours after fully charged at 25ml/h.
Battery recharge time	10h with power on, 3h with power off
Power consumption	25V A
Classification	Class I, type CF
Waterproof Level	IP $\times 3$
Mode of operation	continuous
Dimensions	145 \times 100 \times 120 mm (length \times width \times height)
weight	≤ 1.4 kg
Operating conditions	Temperature:5 $^{\circ}$ ~40 $^{\circ}$,relative humidity:10~95%, Atmosphere pressure:86~106kPa
Transport conditions	Temperature-20 $^{\circ}$ ~+60 $^{\circ}$,relative humidity:10~95%, Atmosphere pressure:50~106kPa
Storage conditions	Temperature-20 $^{\circ}$ ~+45 $^{\circ}$,relative humidity:10~95%, Atmosphere pressure:50~106kPa



Shenzhen Hawk Optical Electronic Instrument Co.,Ltd

Address: 2B, Building No.2, Aerospace Science & Industry factory areas,
Yousong community, Longhua street, Bao'an District, Shenzhen, China.

Postal code: 518109

Tel: 0086-755-8315 1901

Fax:0086-755-8315 1906

Email: szhk@hawkmedical.cn

Web: www.hawkmedical.cn